


VSA 300 S

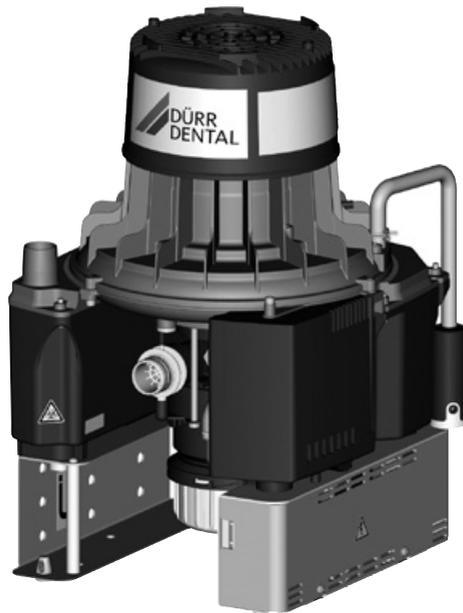

EN

Installation and operating instructions

CE 0297

9000-606-31/30



 **DÜRR
DENTAL**

1905V003

Contents



Important information

1 About this document	3
1.1 Warnings and symbols	3
1.2 Copyright information	3
2 Safety	4
2.1 Intended purpose	4
2.2 Intended use	4
2.3 Improper use	4
2.4 Systems, connection with other devices	4
2.5 General safety information	4
2.6 Specialist personnel	4
2.7 Electrical safety	4
2.8 Only use original parts	5
2.9 Transport	5
2.10 Disposal	5



Product description

3 Overview	6
3.1 Scope of delivery	6
3.2 Accessories	6
3.3 Optional accessories	6
3.4 Consumables	6
3.5 Wear parts and replacement parts	6
4 Technical data	7
4.1 Type plate	11
4.2 Evaluation of conformity	11
4.3 Approvals	11
5 Operation	12
5.1 Tyscor Pulse (optional)	13



Assembly

6 Requirements	14
6.1 Installation/setup room	14
6.2 Setup options	14
6.3 Use of nitrous oxide	14
6.4 Pipe materials	14
6.5 Hose materials	14
6.6 Installation and routeing of hoses and pipes	15
6.7 Information about electrical connections	15
6.8 Information about connecting cables	15
7 System components	16
7.1 Rinsing unit	16
7.2 Spittoon valve	16
7.3 Flow accelerator	16
7.4 Exhaust air filter	16
7.5 Noise reduction	16
8 Installation	17
8.1 Installation and routeing of hoses and pipes	17
8.2 Rinsing unit water connections	19
8.3 Electrical connections	20
8.4 Connections and displays of the control	21
8.5 Display panel connection	22
8.6 Network connection	22
9 Commissioning	23
9.1 Monitoring the unit via the network	23
10 Adjustment options	25
10.1 Setting the afterrun delay time	25
11 Service program	26
12 Description of the service program	27
12.1 Service program ON/OFF	27
12.2 Display test	27
12.3 Sediment level measurement	27
12.4 Motor start - motor braking	27
12.5 Input and output signals	27



Usage

13 Display/handling	28
13.1 Ready for operation	28
13.2 Amalgam collector vessel is 95% full	28
13.3 Amalgam collector vessel is 100% full	28
13.4 Amalgam collector vessel not in position	28
13.5 Motor fault	29
13.6 Brake monitoring	29
14 Monitoring the unit via the network	29
14.1 Monitoring operation	29
14.2 Querying messages	29
14.3 Creating a report	30
15 Disinfection and cleaning	30
15.1 After every treatment	30
15.2 Daily after the end of treatment	31
15.3 Once or twice a week before the midday break	31
16 Replace the amalgam collector vessel	31
16.1 Disposal of amalgam collecting container	32
17 Maintenance	33
17.1 Tests	34



Troubleshooting

18 Tips for operators and service technicians	35
19 Transporting the unit	37

EN



Important information

1 About this document

These installation and operating instructions represent part of the unit.

 If the instructions and information in these installation and operating instructions are not followed, Dürr Dental will not be able to offer any warranty or assume any liability for the safe operation and the safe functioning of the unit.

1.1 Warnings and symbols

Warnings

The warnings in this document are intended to draw your attention to possible injury to persons or damage to machinery.

The following warning symbols are used:

-  General warning symbol
-  Warning – dangerous high voltage
-  Warning – hot surfaces
-  Warning - automatic start-up of the unit
-  Biohazard warning

The warnings are structured as follows:

 **SIGNAL WORD**
Description of the type and source of danger
 Here you will find the possible consequences of ignoring the warning
 > Follow these measures to avoid the danger.

The signal word differentiates between four levels of danger:

- **DANGER**
Immediate danger of severe injury or death
- **WARNING**
Possible danger of severe injury or death
- **CAUTION**
Risk of minor injuries
- **NOTICE**
Risk of extensive material/property damage

Other symbols

These symbols are used in the document and on or in the unit:

-  Note, e.g. specific instructions regarding efficient and cost-effective use of the unit.
-  Refer to Operating Instructions.
-  Wear protective gloves.
-  Wear protective goggles.
-  Disconnect all power from the unit.
-  CE labelling with the number of the notified body
-  Order number
-  Serial number
-  Medical device
-  Manufacturer

1.2 Copyright information

All circuits, processes, names, software programs and units mentioned in this document are protected by copyright.

The Installation and Operating Instructions must not be copied or reprinted, neither in full nor in part, without written authorisation from Dürr Dental.

2 Safety

Dürr Dental has designed and constructed this unit so that when used properly and for the intended purpose it does not pose any danger to people or property. Nevertheless, residual risks can remain. You should therefore observe the following notes.

2.1 Intended purpose

The suction machine provides a suction output to the dental treatment unit and is designed for continuous separation of liquids and air and for separation of amalgam from the entire waste water from dental treatment units.

2.2 Intended use

The suction machine is designed for installation downstream of the hose manifold and spittoon bowl of dental treatment units.

Service, maintenance, recurring tests and cleaning must be performed in accordance with the manufacturer's specifications.

The permissible flow rate must be observed.

The disposable amalgam containers must only be used once.

A rinsing unit is required for surgical procedures and for procedures using prophylaxis powders.

2.3 Improper use

Any use of this appliance / these appliances above and beyond that described in the Installation and Operating Instructions is deemed to be incorrect usage. The manufacturer cannot be held liable for any damage resulting from incorrect usage. The operator will be held liable and bears all risks.

This includes:

- › Use for suction and separation of dust, sludge or plaster.
- › Use in conjunction with flammable or explosive mixtures.
- › Cleaning and disinfection with agents containing sodium hypochlorite or potassium hypochlorite.
- › Assembly in a manner that does not comply with the assembly instructions, in particular installation in rooms containing a potentially explosive atmosphere.
- › Installation without a downward gradient to allow the water to drain off.

2.4 Systems, connection with other devices

Additional devices connected with medical electrical devices must be proven to conform with their corresponding IEC or ISO standards. All configurations must continue to comply with the standard requirements for medical systems (see IEC 60601-1).

Whoever connects additional devices to medical electrical devices automatically becomes the system configurator and is responsible for ensuring that the system corresponds with the standard requirements for systems. Local laws take priority over the requirements outlined above.

2.5 General safety information

- › Always comply with the specifications of all guidelines, laws, and other rules and regulations applicable at the site of operation for the operation of this unit.
- › Check the function and condition of the unit prior to every use.
- › Do not convert or modify the unit.
- › Comply with the specifications of the Installation and Operating Instructions.
- › The Installation and Operating Instructions must be accessible to all operators of the unit at all times.

2.6 Specialist personnel

Operation

Unit operating personnel must ensure safe and correct handling based on their training and knowledge.

- › Instruct or have every user instructed in handling the unit.

Installation and repairs

- › Installation, readjustments, alterations, upgrades and repairs must be carried out by Dürr Dental or by qualified personnel specifically approved and authorized by Dürr Dental.

2.7 Electrical safety

- › Comply with all the relevant electrical safety regulations when working on the unit.
- › Never touch the patient and unshielded plug connections on the device at the same time.
- › Replace any damaged cables or plugs immediately.

Observe the EMC rules concerning medical devices

- › The unit is intended for use in professional healthcare facilities (in accordance with IEC 60601-1-2). If the appliance is operated in another environment, potential effects on electromagnetic compatibility must be taken into account.
- › Do not operate the unit in the vicinity of HF surgical instruments or MRT equipment.
- › Maintain a minimum distance of at least 30 cm between the unit and other electronic devices.
- › Note that cable lengths and cable extensions have effects on electromagnetic compatibility.
- › No maintenance measures are required to maintain the EMV basic safety.

**NOTICE**
Negative effects on the EMC due to non-authorized accessories

- › Use only Dürr Dental parts or accessories specifically approved by Dürr Dental.
- › Using any other accessories may result in increased electromagnetic interference emissions or the unit having reduced electromagnetic immunity, leading to an erroneous operation mode.

**NOTICE**
Erroneous operation mode due to use immediately adjacent to other devices or with other stacked devices

- › Do not stack the unit together with other devices.
- › If this is unavoidable, the unit and other devices should be monitored in order to ensure that they are working correctly.

**NOTICE**
Reduced performance characteristics due to insufficient distance between unit and portable HF communication devices

- › Keep a distance of at least 30 cm between the unit (including parts and cables of the unit) and portable HF communication devices (wireless units) (including their accessories such as antenna cables and external antennas).

2.8 Only use original parts

- › Only use Dürr Dental parts or accessories and special accessories specifically approved by Dürr Dental.
- › Only use only original wear parts and replacement parts.

2.9 Transport

The original packaging provides optimum protection for the unit during transport.

If required, original packaging for the unit can be ordered from Dürr Dental.



Dürr Dental will not accept any responsibility or liability for damage occurring during transport due to the use of incorrect packaging, even where the unit is still under guarantee.

- › Only transport the unit in its original packaging.
- › Keep the packing materials out of the reach of children.

2.10 Disposal

The unit may be contaminated. Instruct the company disposing of the waste to take the relevant safety precautions.

- › Decontaminate potentially contaminated parts before disposing of them.
- › Untamminated parts (e.g. electronics, plastic and metal parts etc.) should be disposed of in accordance with the local waste disposal regulations.
- › If you have any questions about the correct disposal of parts, please contact your dental trade supplier.



3 Overview

3.1 Scope of delivery

The following items are included in the scope of delivery (possible variations due to country-specific requirements and/or import regulations):

VSA 300 S, 230 V AC, 50 Hz 7125-01/002

VSA 300 S, 230 V AC, 50 Hz with in-stalled rinsing unit. 7125-03/002

VSA 300 S, 230 V AC, 60 Hz 7125-04/002

- Set of connection fittings
- Suction hose LW 30, grey
- Hose LW 20
- Hose LW 30, aluminium
- Disposable amalgam container
- Installation and operating instructions
- Operating Handbook
- OroCup

3.2 Accessories

The following items are required for operation of the device, depending on the application:

- Rinsing unit conversion set for VSA 300 S 7100-120-51
- Rinsing unit II 7100-250-50

3.3 Optional accessories

The following optional items can be used with the device:

- Noise reduction hood 7122200000
- Ventilation kit for cabinet installation. 7122-981-51
- Bacteria filter with accessories . . . 7120-143-00

3.4 Consumables

The following materials are consumed during operation of the device and must be ordered separately:

- Disposable amalgam container . . . 7110-033-00
- Orotol plus (2.5 litre bottle) CDS110P6150
- MD 550 spittoon bowl cleaner (750 ml bottle). CCS550C4500
- MD 555 cleaner (2.5 litre bottle) CCS555C6150

3.5 Wear parts and replacement parts

The following working parts need to be changed at regular intervals (refer to the "Maintenance" section):



Information about replacement parts is available from the portal for authorised specialist dealers at: www.duerrdental.net.

4 Technical data

Electrical data		7125-01	7125-04
		7125-03	
Nominal voltage	V	230, 1~	230, 1~
Electrical frequency	Hz	50	60
Nominal current	A	2.9	3.4
Motor protection		Motor winding overheat protector 160 °C (±5 °C)	
Rated power	W	580	770
Type of protection			IP 21
Protection class			I
Protective low voltage	V		24 ~
Output	VA		4

Connections		
Suction connection, DürrConnect special	mm	Ø 30
Exhaust air connection (external)	mm	Ø 30
Drain connection, DürrConnect	mm	Ø 20

Media			
Max. number of operators		1	
Usable volume of collecting container, approx.	ccm	150	
Replacement interval	Months	6 - 9	
Max. unimpeded flow rate	l/min	670	770
Max. suction system pressure	mbar/hPa	-180	-190
Flow rate			
min.	l/min		0.1
max.	l/min		4
Max. suction height	cm		50

General data			
Speed	min ⁻¹	2750	3120
Duty cycle	%		100 (S1)
Dimensions (H x W x D) *	cm	48 x 31 x 31	
Weight, approx.			
without housing	kg		14
with housing	kg		21.5
Noise level ** approx.			
without housing	dB(A)	63	65
with housing	dB(A)	54	56

* Values without accessories and add-on parts

** Noise levels in acc. with EN ISO 1680 "Noise emissions"; measured in a sound-proofed room. The levels are average values with a tolerance of ±3 dB(A). Higher values may be obtained in rooms with reverberating sound characteristics.

Network connection

LAN technology	Ethernet
Standard	IEEE 802.3u
Data rate	Mbit/s 100
Connector	RJ45
Type of connection	Auto MDI-X
Cable type	≥ CAT5

Ambient conditions during storage and transport

Temperature	°C	-10 to +60
Relative humidity	%	< 95

Ambient conditions during operation

Temperature	°C	+10 to +40
Relative humidity	%	< 70

Classification

Medical Devices Directive (93/42/EEC)	Class IIa
---------------------------------------	-----------

Electromagnetic compatibility (EMC)

Interference emission measurements

Interference voltage at the power supply connection CISPR 11:2015/AMD1:2016	Group 1 Class B
Electromagnetic interference radiation CISPR 11:2015/AMD1:2016	Group 1 Class B
Emission of harmonics IEC 61000-3-2:2018	Compliant
Voltage changes, voltage fluctuations and flicker emissions IEC 61000-3-3:2013/AMD1:2017	Compliant

Electromagnetic compatibility (EMC)

Interference immunity measurements cover

Immunity to interference, discharge of static electricity IEC 61000-4-2:2008 ± 8 kV contact ± 2 kV, ± 4 kV, ± 8 kV, ± 15 kV air	Compliant
Immunity to interference, high-frequency electromagnetic fields IEC 61000-4-3:2006+A1:2007+A2:2010 3 V/m 80 MHz - 2.7 GHz 80 % AM at 1 kHz	Compliant
Immunity to interference from power frequency magnetic fields IEC 61000-4-8:2009 30 A/m at 50 Hz	Compliant

EN

Electromagnetic compatibility (EMC)**Interference immunity measurements cover**

Immunity to interference, near fields of wireless HF communication devices

Compliant

IEC 61000-4-3:2006+A1:2007+A2:2010

Immunity to interference table, near fields of wireless HF communication devices

Radio service	Frequency band MHz	Test level V/m
TETRA 400	380 - 390	27
GMRS 460		
FRS 460	430 - 470	28
LTE band 13, 17	704 - 787	9
GSM 800/900		
TETRA 800		
iDEN 820	800 - 960	28
CDMA 850		
LTE band 5		
GSM 1800		
CDMA 1900		
GSM 1900	1700 - 1990	28
DECT		
LTE band 1, 3, 4, 25		
UMTS		
Bluetooth		
WLAN 802.11 b/g/n	2400 - 2570	28
RFID 2450		
LTE band 7		
WLAN 802.11 a/n	5100 - 5800	9

Electromagnetic compatibility (EMC)**Interference immunity measurements supply input**

Immunity to interference, rapid transient bursts – AC voltage grid

IEC 61000-4-4:2012

Compliant

± 2 kV

100 kHz repetition frequency

Immunity to interference, surges

IEC 61000-4-5:2014/AMD1:2017

Compliant

± 0.5 kV, ± 1 kV, L - N

± 0.5 kV, ± 1 kV, ± 2 kV, L/N - PE

Electromagnetic compatibility (EMC)

Interference immunity measurements supply input

Immunity to interference, line-conducted disturbances induced by high-frequency fields – AC voltage grid

IEC 61000-4-6:2013

3 V

0.15 - 80 MHz

6 V

ISM frequency bands

6.765 - 6.795 MHz

13.553 - 13.567 MHz

26.957 - 27.283 MHz

40.66 - 40.70 MHz

80 % AM at 1 kHz

Compliant

Immunity to interference due to voltage dips, short interruptions and voltage variations

IEC 61000-4-11:2004/AMD1:2017

Compliant

Electromagnetic compatibility (EMC)

Interference immunity measurements SIP/SOP

Immunity to interference, rapid transient bursts – I/O, SIP/SOP ports

IEC 61000-4-4:2012

± 1 kV

100 kHz repetition frequency

Compliant

Immunity to interference, line-conducted disturbances induced by high-frequency fields – SIP/SOP ports

IEC 61000-4-6:2013

3 V

0.15 - 80 MHz

6 V

ISM frequency bands

6.765 - 6.795 MHz

13.553 - 13.567 MHz

26.957 - 27.283 MHz


40.66 - 40.70 MHz

80 % AM at 1 kHz

Compliant

4.1 Type plate

The type plate is located on the noise reduction hood.

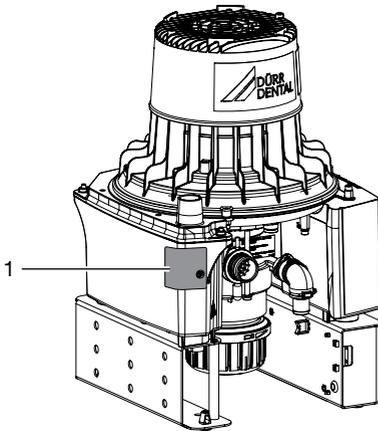

1 Type plate

4.2 Evaluation of conformity

This device has been subjected to conformity acceptance testing in accordance with the current relevant European Union guidelines. This equipment conforms to all relevant requirements.

4.3 Approvals


**Centre of Competence in Civil Engineering,
Berlin**

Test number Z-64.1-15

Separation method compliant with standard

ISO 11143 Type 1

5 Operation

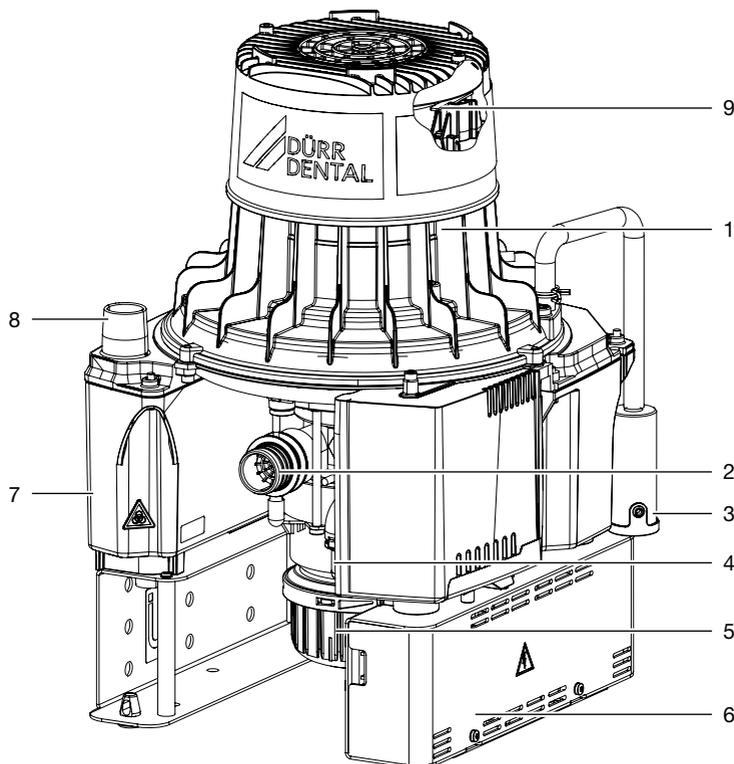

- 1 Motor
- 2 Inlet connection with protective coarse filter
- 3 Auxiliary air nozzle
- 4 Waste water connection
- 5 Amalgam collecting container
- 6 Control electronics
- 7 Exhaust air muffler
- 8 Exhaust air connection
- 9 Rotational speed sensor

The mixture of liquids, solid particles and air drawn in passes through the inlet connection and into the suction unit. The coarse filter holds back the solid particles.

Inside the separation unit, the aspirated fluids and solid particles pass through a two-stage separation system and are separated from the suction air. This separation system consists of a cyclone separator and a separation turbine.

The aspirated mixture flows into the cyclonic separator, where it is set into a spiral motion. In this first stage, the resulting centrifugal forces force the fluid constituents and any remaining solid particles against the outside wall of the separation chamber of the cyclone separator. This initially only effects a "coarse separation" of the fluid.

In the subsequent second stage, the separation turbine effects "fine separation" of the remaining liquid from the air flow that has carried it to this point.

The fluid and solids accreting in the separation chamber are continuously fed to the amalgam centrifuge, where the amalgam particles are removed. The fluid extracted via the centrifuge is fed through the waste water valve and the outlet connection into the central waste-water system.

Under the centrifuge there is an interchangeable collecting container into which the separated amalgam particles fall once the motor is switched off.

A sensor checks the fill level in the collecting container, and when it is full, an LED on the display panel indicate that the collecting container needs to be replaced. Depending on the type of work carried out and the amount of amalgam arising, the collector vessel should be changed approx. every 6-9 months. A secure twist cap makes the replacement and closing of the collector vessel easier.

A pump connected to the centrifuge keeps the fluid level constant in the collecting container. This prevents accidental overflow when replacing the collecting container.

The air separated from the liquid is sucked off by the vacuum pressure generated by the turbine wheel. The air is then blown through the noise reduction hood and over the exhaust air connection and out of the machine.

The turbine wheel, separation turbine and amalgam centrifuge are driven by the motor.

An auxiliary air nozzle is connected to the turbine housing. The auxiliary air nozzle limits the vacuum in the system. In certain working situations it also sucks additional cooling air into the machine.

5.1 Tyscor Pulse (optional)

The software is connected via the network to the devices from Dürr Dental and displays the current status as well as messages and errors.

All messages are logged and can be printed or sent.

The *cockpit* shows the devices with the current characteristic data and provides a quick overview of the functional status of the devices.



6 Requirements

The unit can be installed on the same level as the surgery room or in a floor below.



Further information can be found in our suction planning information leaflet. Order number 9000-617-03/..

6.1 Installation/setup room

The room chosen for set up must fulfil the following requirements:

- Closed, dry, well-ventilated room
- Should not be a room made for another purpose (e. g. boiler room or wet cell)



Ambient and environmental conditions must be taken into account. Do not operate the unit in damp or wet conditions.

- When installing in a cabinet the inlet and outlet ventilation slots must be present; minimum free cross-section at least 120 cm².
- Forced ventilation (fan) must be provided if there is a risk that the recommended room air temperature could be exceeded. The air flow performance must be at least 2 m³/min.
- Do not cover cooling slots or openings with housing installations; ensure sufficient clearance to the openings to permit sufficient cooling.

6.2 Setup options

The following options for setting up the unit are available:

- Wall installation using a Dürr Dental wall mounting
- In a ventilated cabinet
- In a Dürr Dental noise reducing housing



When installed in a basement or a similar room, the unit must be placed on a base or be fixed to the wall at a minimum height of 30 cm above the floor.

6.3 Use of nitrous oxide

This unit is technically suitable for the aspiration of nitrous oxide (laughing gas). However, when assembling a system for aspiration of nitrous oxide, it is important to ensure that the other components in the system are also suitable for this purpose. Those responsible for setting up the system must assess this and approve and release the system for the aspiration of nitrous oxide.



Operation with nitrous oxide is only permitted if the exhaust air is transported from the unit to the outside of the building.

6.4 Pipe materials

Only use vacuum-sealed HT-waste pipes manufactured from the following materials:

- Polypropylene (PP),
- Chlorinated polyvinyl chloride (PVC-C),
- Unplasticized polyvinyl chloride (PVC-U),
- Polyethylene (PE).

The following materials must not be used:

- Acrylonitrile-butadiene-styrene (ABS),
- Styrene copolymer blends (e.g. SAN + PVC).

6.5 Hose materials

For waste connections and suction lines only use the following hose types:

- Flexible spiral hoses made of PVC with integrated spiral or equivalent hoses
- Hoses that are resistant to dental disinfectants and chemicals



Plastic hoses will display signs of ageing over time. Therefore, they should be inspected regularly and replaced as necessary.

The following types of hoses must not be used:

- Rubber hoses
- Hoses made completely of PVC
- Hoses that are not sufficiently flexible

6.6 Installation and routing of hoses and pipes

- › Execute the on-site pipe installation in accordance with the applicable local regulations and standards.
- › Lay the hose installation of the drains to or from the unit at a sufficient incline.



If incorrectly laid, the hoses can become blocked with sedimentation.

6.7 Information about electrical connections

- › Ensure that electrical connections to the mains power supply are carried out in accordance with current valid national and local regulations and standards governing the installation of low voltage units in medical facilities.
- › Install an all-pole disconnect switch with a contact opening width of at least 3 mm in the electrical connection to the mains power supply.
- › Observe the current consumption of the devices that are to be connected.

Electrical fusing

LS switch 16 A, characteristic B, C and D in accordance with 60898.

6.8 Information about connecting cables

The diameter of the connections depends on the current consumption, length of line and the ambient temperature of the unit. Information concerning the current consumption can be found in the Technical Data supplied with the particular unit to be connected.

The following table lists the minimum diameters of the connections in relation to the current consumption:

Current consumption of unit [A]	Cross-section [mm ²]
> 10 and < 16	1.5
> 16 and < 25	2.5
> 25 and < 32	4
> 32 and < 40	6
> 40 and < 50	10
> 50 and < 63	16

Mains supply cable

Installation type	Line layout (minimum requirements)
Fixed installation	– Plastic sheathed cable (e.g. type NYM-J)
Flexible	– PVC flexible line (e.g. H05 VV-F) or – Rubber connection (e.g. H05 RN-F or H05 RR-F)

Display panel

Installation type	Line layout (minimum requirements)
Fixed installation	– CAT5.e network cable
Flexible	– ISDN standard cable with connectors or – Network patch cable

Control cable

24 V protective low voltage for:

- Hose manifold
- Place selection valve
- Spittoon valve

Installation type	Line layout (minimum requirements)
Fixed installation	– Shielded sheathed cable (e.g. (N)YM (St)-J)
Flexible	– PVC data cable with shielded cable sheathing, as used for telecommunications and IT processing systems (e.g. type LiYCY) or – Lightweight PVC control cable with shielded cable sheathing

7 System components

The system components listed below are required or recommended for various procedures or for installation.

7.1 Rinsing unit

In the absence of a spittoon or rinsing unit in the treatment unit, it is absolutely essential that a rinsing unit be installed in the VSA 300 S. In addition, for surgical procedures and for procedures using airflow a rinsing unit must always be installed in the treatment unit to supply a small amount of water to the system during aspiration. Any secretions present will thus be diluted and can be transported away more easily.

For further information, refer to the rinsing unit installation and operating instructions

7.2 Spittoon valve

To ensure that the waste water from the spittoon bowl is directed into the suction pipe, a spittoon valve needs to be installed between the drain of the spittoon bowl and the suction pipe. The spittoon valve must switch on the suction machine if required.

7.3 Flow accelerator

In order to keep the suction system free of deposits, a flow accelerator can be fitted in conjunction with a spittoon valve. When using a bowl rinse system, water will collect before the flow accelerator. The next time suction takes place using the large cannula, the collected fluid is transported in surges and at high speed to the suction system. This ensures automatic cleaning of the suction pipes.

7.4 Exhaust air filter

For hygienic reasons, we recommend the installation of a bacteria filter in the exhaust air line.

If the unit is installed in the surgery and the exhaust air cannot be discharged to the outdoors, it is essential to install a bacteria filter.

Depending on the type and condition of the bacteria filter, it will need to be replaced every 1-2 years at the latest.



The separation integrated in the system does not retain bacteria; this is why we recommend installing a suitable filter in the exhaust air system.

7.5 Noise reduction


If the noise level from the exhaust air vent or the flow noise generated is too high, noise reduction can be installed in the exhaust air line.

8 Installation

 The actual connection can vary depending on the chosen installation option. The connection shown is only an example.

8.1 Installation and routing of hoses and pipes

- › Establish connections between the pipe system and the unit using the flexible hoses supplied. This will prevent vibrations from being transmitted to the pipe system.
- › The connection between the pipe line and unit suction connection should be kept as short as possible and straight, without bends.
- › Install the drain hoses with a downward gradient so that the waste water can drain off.
- › Waste water connections must be implemented in accordance with applicable local and national regulations.

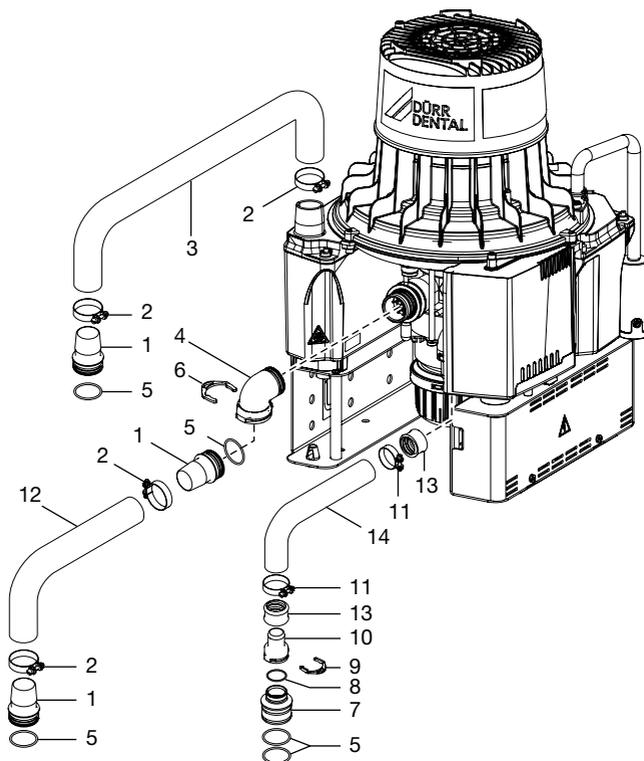

- 1 Hose connector \varnothing 30 mm
- 2 Hose clip 25-40 mm
- 3 Waste air pipe (aluminium) \varnothing 30 mm inside
- 4 Elbow DN 30
- 5 O-ring \varnothing 30x2 mm
- 6 Securing ring
- 7 Connector \varnothing 36 mm (external)
- 8 O-ring \varnothing 20x2 mm
- 9 Securing ring
- 10 Hose connector socket \varnothing 20 mm




- 11 Hose clip \varnothing 28 mm
- 12 Suction hose \varnothing 30 mm (internal)
- 13 Hose sleeve
- 14 Waste water hose \varnothing 20 mm (internal)


8.2 Rinsing unit water connections

 Check the water pressure for the rinsing unit. The water pressure should be between 2 and 4 bar.

- › Screw the Tecalan hose with sleeve piece, double-tapered ring and locking nut onto the rinsing unit.

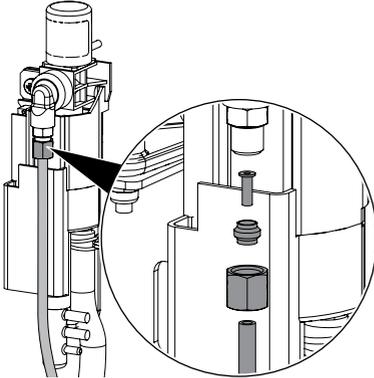

- › Apply the T-piece for Tecalan water hose with $\varnothing 4$ mm or $\varnothing 6$ mm in the water supply.
- › Apply the Tecalan hose with sleeve piece, double-tapered ring and locking nut to the T-piece.

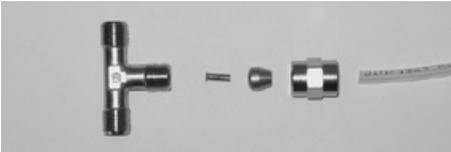


- › Alternatively, apply the Tecalan hose with adapter piece, seal, R3/4" screw connection, sleeve piece, double-tapered ring and locking nut onto a water tap.




8.3 Electrical connections

- › Connect the control line.
- › Connection the display panel.
- › Connect the network cable (optional when using monitoring software).
- › Establish the electrical connection to the supply network (230 V).

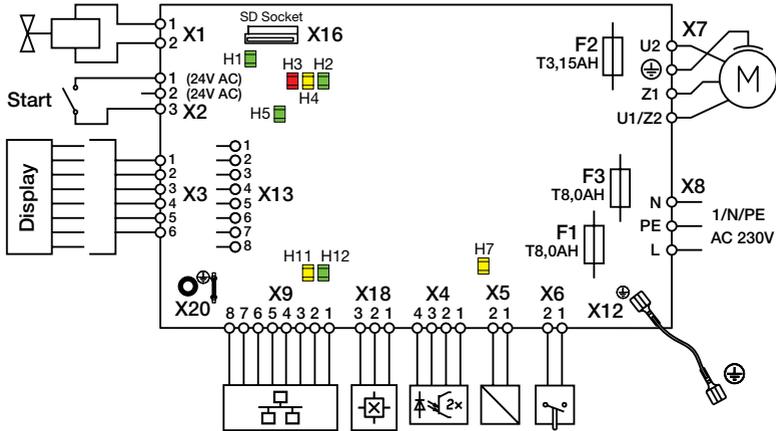

- X1 Voltage supply for the rinsing unit
- X2 24V output voltage and switching contact to suction unit in the treatment unit
- X3 Display panel
- X4 Sediment sensor light barriers
- X5 Sediment sensor lifting magnet
- X6 Collecting container safety switch
- X7 Motor connection
- X8 Mains connection
- X9 Network connection
- X12 Ground contact to the unit housing
- X13 Display panel
- X16 SD card holder (for Micro SD)
- X18 Connection of the hall sensor for monitoring the speed
- X20 Ground contact to the unit housing
- F1 Main fuse
- F2 Brake fuse
- F3 Main fuse
- H1 Rinsing unit
- H2 Display green (as with display panel)
- H3 Display red (as with display panel)
- H4 Display yellow (as with display panel)
- H5 Switching contact control signal suction unit in the treatment unit
- H7 Sediment sensor lifting magnet
- H11 Network
- H12 Network

8.4 Connections and displays of the control

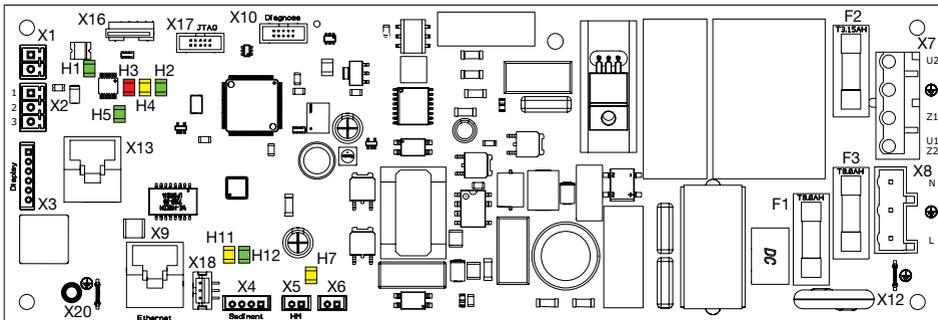

- X1 Voltage supply for the rinsing unit
- X2 24V output voltage and switching contact to suction unit in the treatment unit
- X3 Display panel
- X4 Sediment sensor light barriers
- X5 Sediment sensor lifting magnet
- X6 Collecting container safety switch
- X7 Motor connection
- X8 Mains connection
- X9 Network connection
- X10 Diagnosis
- X12 Ground contact to the unit housing
- X13 Display panel
- X16 SD card holder (for Micro SD)
- X17 JTAG programming interface
- X18 Connection of the hall sensor for monitoring the speed
- X20 Ground contact to the unit housing
- F1 Main fuse
- F2 Brake fuse
- F3 Main fuse
- H1 Rinsing unit
- H2 Display green (as with display panel)
- H3 Display red (as with display panel)
- H4 Display yellow (as with display panel)
- H5 Switching contact control signal suction unit in the treatment unit
- H7 Sediment sensor lifting magnet
- H11 Network
- H12 Network


8.5 Display panel connection




There must be a direct line connecting the network socket on the unit and the network socket on the display panel. Do not toggle network units (e. g. switch or router).

Observe the resistance of the network cable between the network sockets. The maximum length should not exceed 50 m.

- Connect the network cable in the network socket to the VSA 300 S in the network bushing (X13).
- Connect the display panel and network socket using the ISDN cable supplied.

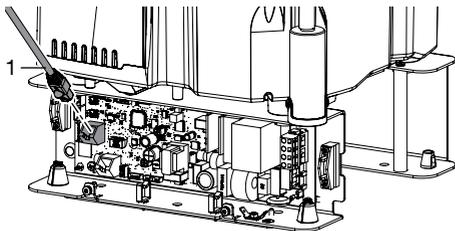

1 Network cable

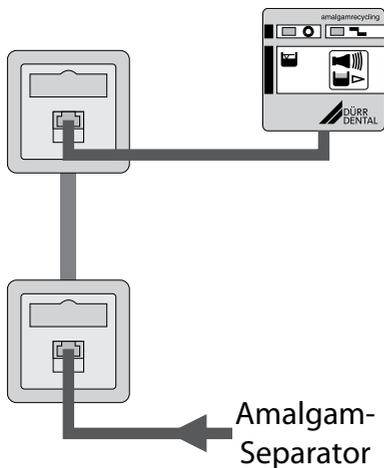


8.6 Network connection

Purpose of the network connection

The network connection is used to exchange information or control signals between the unit and a software installed on a computer, in order to, e. g.:

- Display parameters
- Select operating modes
- Indicate messages and error situations
- Change unit settings
- Activate test functions
- Transmit data for archiving
- Provide documents concerning the units
- Plug in the network cable at the network connection on the unit (optional when using monitoring software).
- Plug in the network cable at the network socket.

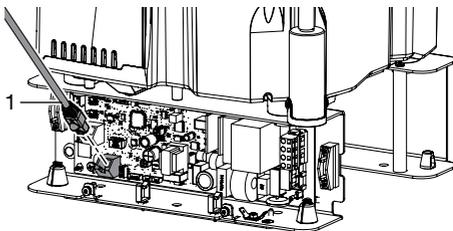

1 Network cable

9 Commissioning



NOTICE

Interference caused by larger particles such as pieces of tooth or fillings

- › Do not operate the unit without a coarse filter

- › Check that the coarse filters are installed in the suction system (e.g. in the spittoon).
- › Turn on the unit power switch or the main surgery switch.
- › Carry out a function check of the device.
- › Check all connections for leak tightness.
- › Carry out an electrical safety check in accordance with applicable regulations (e.g. regulations concerning set up, operation and application of medical devices) and record the results as appropriate, e.g. in the technical log book.



In many countries technical medical products and electrical devices are subject to regular checks at set intervals. The owner must be instructed accordingly.

9.1 Monitoring the unit via the network

Combining devices safely

- The overall safety of the unit and its main performance features are independent of the network. The device is designed for operation independent of a network. However, some of the functions are not available in this case.
- Incorrect manual configuration can lead to significant network problems. The expert knowledge of a network administrator is required for configuration.
- The data connection utilises part of the bandwidth of the network. Interactions with other medical devices cannot be completely ruled out. Apply the IEC 80001-1 standard for risk assessment.
- The device is not suitable for direct connection to the public internet.
- › When connecting the unit to other devices, such as a PC system, comply with the requirements set out in section 16 of IEC 60601-1 (EN 60601-1).

- › When setting up the PC system in the vicinity of the patients:
 - Only connect components (e.g. computer, monitor, printer) that comply with the standard IEC 60601-1 (EN 60601-1).
- › When setting up the PC system outside of the vicinity of the patients:
 - Connect components (e.g. computer, monitor, printer) that comply at least with the standard IEC 60950-1 (EN 60950-1) at least.

Network configuration

Various options are available for network configuration:

- Automatic configuration via DHCP (recommended).
- Automatic configuration via Auto-IP for direct connection of unit and computer.
- Manual configuration.
- › Configure the network settings of the unit using the software or, if available, the touch screen.
- › Check the firewall and release the ports, if applicable.

Network protocols and ports


Port	Purpose	Service
45123 UDP, 45124 UDP	Unit recognition and configuration	
1900 UDP	Service indicator	SSDP / UPnP
502 TCP	Unit data	
514 ¹⁾ UDP	Event protocol data	Syslog
22 TCP, 23 TCP	Diagnosis	Telnet, SSH
123 UDP	Time	NTP

- ¹⁾ The port can vary depending on the configuration.

The following requirements must be met in order to monitor the unit on the computer:

- Unit connected to the network
- Current monitoring software installed on the computer

Add device

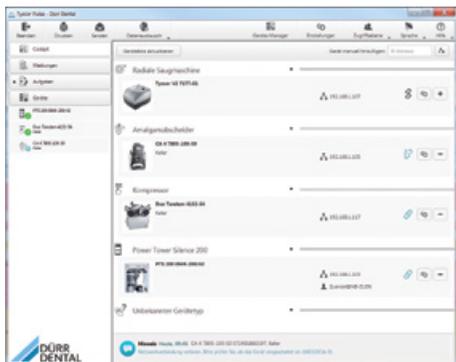

Connection icons:



The unit is present in the network and connected to the software.



The unit is present in the network but not connected to the software.



The network connection between the software and the unit has been interrupted, e.g. the device is switched off.

Requirements:


- Unit is switched on and connected to the network
- Administrator or service technician access level selected in the software
- > Working in the menu bar, click on  **Device Manager**.

The list of units appears. A symbol displays the connection status to the software:

The new unit that is not yet connected, is displayed with the connection status .

- > Select the unit and click on .
- The unit appears in the side bar.

Adding the device in the cockpit

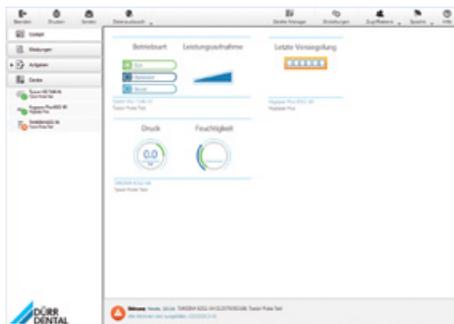

All devices that are connected to the software can be added to the cockpit. When the unit is first connected to the software, the unit is automatically added to the cockpit.


Requirements:

- Administrator or Service Technician access level selected.
- > Click on the device in the device list with the left mouse button and keep the mouse button pressed.
- > With the mouse key pressed, drag the unit onto the cockpit.
- > Release the mouse key.

The block with the current characteristic data and the name of the device appear in the cockpit.

- > To change the position of the device block, click on the block and, with the mouse key pressed, drag it to the required location.

Manually starting the device

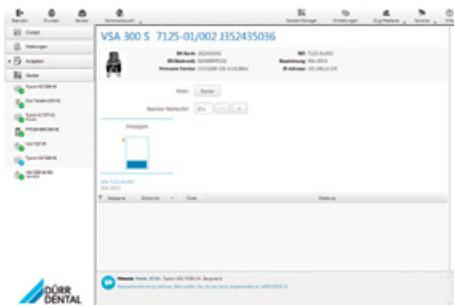

Manually starting the device for testing.

Requirements:

- Service technician access level selected.
- › Select the device in the device list.
- › Click on the Start button with the left mouse key; on some devices you will need to keep it pressed.

10 Adjustment options

10.1 Setting the afterrun delay time



The delay time can be set via Tyscor Pulse.

Requirements:

- Administrator or Service Technician access level selected.
- › Select the device in the list of devices.
The block with the current button parameters and the name of the device appear in the Contents area.
- › Use the "+" and "-" buttons to adjust the time in increments of 1 s.

11 Service program

START

1

2

3

4

STOP

12 Description of the service program



Wear protective equipment to avoid any risk of infection (e.g. liquid-tight protective gloves, protective goggles, face mask).

The various unit functions can be checked with the aid of the service program.

The individual program steps are:

- Display test
- Sediment level measurement
- Motor start and motor braking with rpm check
- Input and output signals

Function of the service key:

- By pressing the service key twice the next individual program step is called up.
- By pressing the service key once that program step is repeated.

A press of the service key is confirmed by an audible signal.

12.1 Service program ON/OFF

On

- Press the service key and switch on the voltage supply to the unit.
- As soon as a signal melody can be heard, release the service key.

The green, yellow and orange LEDs on the display panel light up (display test) and the service program is activated.

Off

Switch off the main supply to the unit.

12.2 Display test

The display test is activated as soon as the service program is started.

The LEDs on the display panel are checked. All three LEDs must come on. There is also an audible signal, which can be switched off by pressing the service button.

12.3 Sediment level measurement



While the service program is activated, the safety check for the collector vessel is deactivated.

The sediment level measurement can be used to check the function of the sediment sensor and the function of the LEDs.

Every time the service key is pressed, the sediment level is checked. If a test container is used, the different 95% and 100% filling level on the display panel can be revealed.

12.4 Motor start - motor braking

The drive motor starts up and is automatically braked after the delay time. If the service key is pressed before the end of the delay time, the motor will immediately be braked.

This procedure can be repeated by pressing the service key 1x again.


As a result of the rpm monitoring, the LED will go from orange to green on start-up and from green to orange during braking.

12.5 Input and output signals

- After activating the program point, the yellow LED on the display panel flashes.
- A cycled DC voltage (c. 22-30 V) can be measured on the rinsing unit connection (X1).
- Opening the collecting container causes the orange display to illuminate on the display panel.
- If a start signal is applied to socket X2 (lift out the suction hose on the hose manifold) the green LED illuminates on the display panel.



13 Display/handling

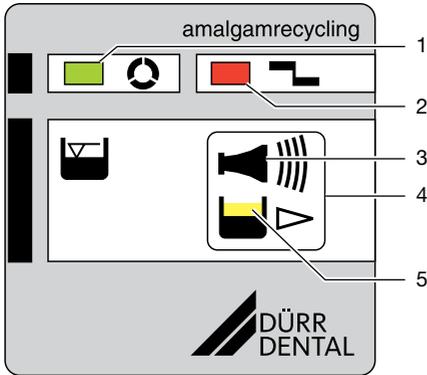

- 1 GREEN LED
- 2 RED display
- 3 Audible signal/melody
- 4 Reset/service key
- 5 YELLOW LED

13.1 Ready for operation

- GREEN LED lights up

13.2 Amalgam collector vessel is 95% full

- Yellow LED lights up
- GREEN LED lights up

- Audible signal melody sounds

- At a fill level of 95%, the signal melody can be switched off by pressing the reset button. The device is then ready for operation again.
- The yellow LED comes on as a reminder that the amalgam collecting container is due to be changed. The level display is repeated every time the unit is switched on at the main power switch.

We recommend changing the amalgam collecting container when it reaches 95% full.

13.3 Amalgam collector vessel is 100% full

- Yellow LED lights up
- Red display flashes
- Audible signal melody sounds

- At a fill level of 100% the signal melody can no longer be switched off by pressing the reset button.

- The collecting container needs to be replaced.

Wear protective equipment to avoid any risk of infection (e.g. liquid-tight protective gloves, protective goggles, face mask).

- The separator will not be ready for operation again until the amalgam collecting container has been replaced

13.4 Amalgam collector vessel not in position

- Red display flashes
- Audible signal

- Press the reset button briefly to switch off the audible signal.
- Switch off the device.
- Insert the collecting container.
- Switch on the unit.
- Green LED lights up – "Ready for operation"

If this error message occurs when the collecting container is correctly inserted, this indicates that there is a technical defect – inform your Service Technician.

13.5 Motor fault

 Red display and

 green LED flash alternately

 Audible signal

 Occurs during the start-up of the amalgam separator.

- Press the reset button briefly to switch off the audible signal.
- If the reset button is pressed for longer than 2 seconds the unit can be restarted.

 If this problem happens again on the same day, the amalgam separator will no longer be operational - notify the service technician.

13.6 Brake monitoring

 Red display and

 green LED flash alternately

 Audible signal

 Occurs upon braking action of amalgam separator.

- Press the reset button briefly to switch off the audible signal.
- The amalgam separator is still operational.


 If this problem occurs on several consecutive days, the braking must be checked by a service technician.

14 Monitoring the unit via the network

 As the monitoring system of the device, the software must deliver acoustic signals. Audio output on the computer must be activated.

14.1 Monitoring operation

The device must have been added to the cockpit for the graphical device block to be shown in the cockpit.

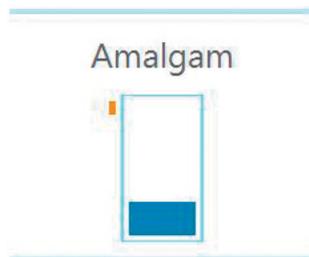

The following is shown in the unit block of the amalgam separator:

- Fill levels in the collecting container

14.2 Querying messages

-  Trouble-free operation
-  Fault
Operation of the device interrupted
-  Warning
Operation of the device restricted
-  Note
Important information about the device
-  Information
-  Establishing a connection to the device
-  Connection to the device interrupted

If a message occurs for an device, the symbol next to the device in the side bar changes. The message appears in the cockpit and in the device details.

If several messages occur, the symbol of the highest message level in each case is displayed.

 As soon as a message concerning a device occurs, the symbol in the task bar (or Mac OS menu bar) also changes to the relevant message symbol. If required by the message an acoustic signal also sounds.

› To query the message details, switch to the cockpit or to the device.

14.3 Creating a report

You can print out a current report  or sent it via e-mail .

The report contains all messages and a screenshot of the view that is displayed when the report is created.

15 Disinfection and cleaning



NOTICE

Device malfunctions or damage due to use of incorrect media

Guarantee claims may become invalid as a result.

- › Do not use any foaming agents such as household cleaning agents or instrument disinfectants.
- › Do not use abrasive cleaners.
- › Do not use agents containing chlorine.
- › Do not use any solvents like acetone.

Dürr Dental recommends

- For disinfection and cleaning:
Orotol plus or Orotol ultra
- For cleaning:
MD 555 cleaner

Only these products have been tested by Dürr Dental.

When using prophylaxis powders, Dürr Dental recommends the water-soluble Lunos prophylaxis powders in order to protect the Dürr Dental suction systems.

15.1 After every treatment

- › Aspirate a glass of cold water through the large and the small suction hoses. Do this even if only the small suction hose was actually used during treatment.



Suction through the large suction hose causes a large amount of air to be drawn up, thereby considerably increasing the cleaning effect.

15.2 Daily after the end of treatment

 After higher workloads before the midday break and in the evening

The following are required for disinfection/cleaning:

- Non-foaming disinfectant/cleaning agent that is compatible with the materials.
- Unit care system, e.g. OroCup
- › To pre-clean, suck up 2 litres of water with the care system.
- › Aspirate the disinfection/cleaning agent with the care system.

15.3 Once or twice a week before the midday break

 Under harsher conditions (e.g. hard water or frequent use of prophylaxis powders) 1x daily before the midday break

The following are required for cleaning:

- Special non-foaming detergent for suction units that is compatible with the materials.
- Unit care system, e.g. OroCup
- › To pre-clean, suck up 2 litres of water with the care system.
- › Aspirate the cleaning agent with the care system.
- › Rinse with ca. 2 l water after the application time.

16 Replace the amalgam collector vessel



WARNING

Risk of contamination if the amalgam collector vessel is reused since the collector vessel is not water-tight.

- › Do not use the collector vessel more than once (disposable item).



Wear protective equipment to avoid any risk of infection (e.g. liquid-tight protective gloves, protective goggles, face mask).



We strongly recommend that the amalgam collecting container should only be changed in the morning before the start of work. This will prevent fluid from dripping out of the drum while it is being changed.

- › Disconnect all power from the unit.
- › Remove the full amalgam collecting container and from the device.
- › Pour disinfectant for suction units (e.g. Orotol plus, 30 ml) into the full amalgam collecting container.
- › Close and secure the full amalgam collecting container using the cap. Observe the markings on the cap and on the collecting container.
- › Place the securely closed amalgam collecting container into its original packaging and seal.
- › Insert a new amalgam collecting container in the unit and lock it in position with the vessel lift.



Only use original amalgam collecting container.

- › Switch on the power supply. The unit is ready for operation again.

16.1 Disposal of amalgam collecting container



The contents of the amalgam collecting container are contaminated with heavy metals and must not be disposed of as household waste or the environment.

- Collection and waste disposal by a waste management company specialised in surgery waste.
- Collection and waste disposal by an approved waste management company.

EN

17 Maintenance



All maintenance work must be performed by a qualified expert or by one of our Service Technicians.



Prior to working on the unit or in case of danger, disconnect it from the mains.



WARNING

Infection due to contaminated unit

- › Clean and disinfect the suction before working on the unit.
- › Wear protective equipment when working (e. g. impermeable gloves, protective goggles and mouth and nose protection).

EN

Maintenance interval	Maintenance work
Dependent upon the level of usage of the device	<ul style="list-style-type: none"> › Replace the amalgam collecting container when a fill level of 95% or 100% is indicated on the display panel <p>Notes concerning prophy powders:</p> <p>The amalgam separator is not functionally affected by conventional prophy powders. Under certain circumstances however, increased soiling of lines and hoses and a more frequent changing of the amalgam collecting container can be expected.</p>
Annually	<ul style="list-style-type: none"> › Cleaning of the suction unit in accordance with the operating instructions. › Check the inlet and outlet hoses for signs of deposits/blockage or cracks and replace where necessary. › Check the outflow valve and replace if necessary. › Replace the exhaust air filter (depending on the installation conditions).

17.1 Tests



WARNING

Infection due to contaminated unit

- › Clean and disinfect the suction before working on the unit.
- › Wear protective equipment when working (e. g. impermeable gloves, protective goggles and mouth and nose protection).



In some countries the owner is required to keep an operating handbook. This operating handbook must document all maintenance work, service work, checks and amalgam disposal.

Annual inspection

This inspection should only be carried out by suitably trained staff.

Work steps to be performed:

- › General functional check (e.g. aspiration, spittoon inlet)
- › During the sediment fill level measurement, visually inspect the operability of the sediment sensor.
- › Service program

Tyscor Pulse (optional)

This test should be performed as an additional test if the device is monitored with Tyscor Pulse.

Requirements for the test:

- Device connected to the network.
- Tyscor Pulse has been started.

Work steps to be performed:

- › Check whether any messages are displayed on the PC monitor.
- › Check the acoustic signal.

Inspection of the general operating condition every 5 years

This inspection must be carried out every 5 years (in accordance with the German Waste Water Regulations, Annex 50, Dental Treatment) by an inspector in accordance with national regulations.

For inspection, the following are required:

- Test vessel
- Measuring beaker

Work steps to be performed:

- › Remove the collector vessel. Here, the red LED on the display panel must flash and an audible signal must sound.
- › Insert the test collector vessel.
- › Press the service key on the display panel.
- › Suck up c. 1 L water.
- › Once the device has switched off, remove the test vessel and measure the remaining amount of water.

The unit is working correctly if:

- there is at minimum content of 70 ml in the test vessel.

If there is less fluid, clean the centrifuge drum or check the operation of the unit.



18 Tips for operators and service technicians



Any repairs exceeding routine maintenance may only be carried out by qualified personnel or our service.



WARNING

Infection due to contaminated unit

- › Clean and disinfect the suction before working on the unit.
- › Wear protective equipment when working (e. g. impermeable gloves, protective goggles and mouth and nose protection).



Prior to working on the unit or in case of danger, disconnect it from the mains.

Fault	Probable cause	Solution	
Device does not start	No mains voltage	› Check the mains supply voltage. * › Check the fuses and replace if necessary. *	
	Undervoltage	› Measure the supply voltage; call an electrician if necessary. *	
	Control electronics defective	› Replace the electronics. *	
Motor fault	Rotational speed sensor not working	› Check that the Hall sensor is correctly seated. * › Check the plug connections of the sensor cable. * › Check the magnets in the fan wheel. *	
	The unit generates unusual noises	Solid particles in the turbine chamber	› Disassemble the unit and clean the turbine and housing. *
	Water leaking from the exhaust air connection	Membrane valve blocked	› Check the membrane valve at the waste water connection and if necessary clean or replace. *
Foam in turbine due to use of incorrect cleaning and disinfectant agents		› Use non-foaming cleaning and disinfectant agents.	
Build-up of condensate in the exhaust air line		› Check the pipe system; avoid over-cooling. *	
Waste water line/siphon trap blocked		› Clean the waste water line/siphon trap. *	

Fault	Probable cause	Solution
Suction performance too low	Coarse filter blocked	› Clean the coarse filter at the intake connection.
	Leak in the suction line	› Check and if necessary establish leak-tightness of suction system and connections. *
	Mechanical sluggishness of turbine caused by soiling	› Disassemble the unit and clean the turbine. *

* Only to be done by service technicians.

19 Transporting the unit



WARNING


Infection due to contaminated unit

- › Disinfect the unit before transport.
- › Close all media connections.



Wear protective equipment to avoid any risk of infection (e.g. liquid-tight protective gloves, protective goggles, face mask).

- › Before disassembly, clean and disinfect the suction unit and the unit using a suitable disinfectant approved by Dürr Dental.
- › Disinfect a defective unit using a suitable surface disinfection agent.
- › Seal all connections with sealing caps.
- › Pack the unit securely in preparation for transport.

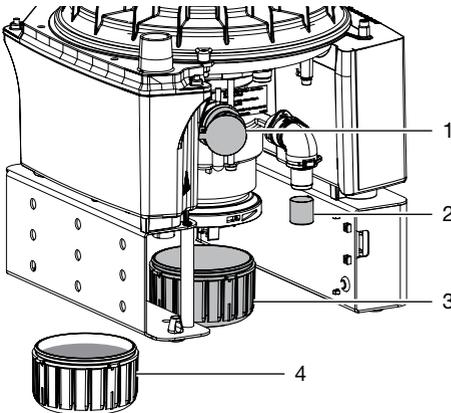

- 1 Suction connection sealing caps
- 2 Water outflow sealing cap
- 3 Amalgam collecting container EMPTY
- 4 Amalgam collecting container

EN



Hersteller/Manufacturer:

DÜRR DENTAL SE
Höpfigheimer Str. 17
74321 Bietigheim-Bissingen
Germany
Fon: +49 7142 705-0
www.duerrdental.com
info@duerrdental.com

